



MedReady[™]

Available with cellular compliance monitoring



medreadyinc.net

Automated Medication Dispenser

MedReady, Inc.

2103 Border Ave. Torrance, CA 90501 310-328-7557 Sales@MedReady.net One Less Argument • More Independence
Feel Better • Great Price

- Timed-alarm system provides next dosage automatically
- Holds up to 7 days of 4 times a day dosing
- MedReady's functions and capabilities provide the best value in the medication dispensing marketplace
- The cellular dispenser has the capability to communicate compliance events via cell towers. No phone line or internet connection is necessary at the dispenser location.
- Check compliance history, change dispensing times, and update all device settings via the www.medready.net website.



FEATURES

Large 1/2" L.E.D. displays the time and aids setting alarms.

- Medication regimen can dose up to 4 times per day:
 - ◆ "Time for medicine" alarm beeps 5 seconds on, 1 second off continuously until the medication door is slid all the way open exposing the current dose of medication.
 - Programmable alarm length sounds for as short as 3 minutes and as long as 250 minutes in 10 minute increments.
 - Powerful, 90dB, 400Hz. audible alarm is designed for maximum effectiveness.
 - Audible and visual alarms can be independently turned on or off.
- MedReady holds 28 doses of dry medication such as pills or capsules.
 - ♦ 4 Doses a day for 7 days
 - ♦ 3 Doses a day for 9 days
 - ◆ 2 Doses a day for 14 days
 - ◆ 1 Dose a day for 28 days
- Pill compartments are "ramped" for easy removal and can fit up to 9 regular man.

MedReady[™]

Automated Medication Dispenser

Available with cellular compliance monitoring

- Only the most recent medication dose is available. Pill tray rotates to the next dose when the next audible alarm sounds. If a patient misses a dose, the pill dispenser denies access to the missed dose.
- Alarm in pill dispenser is as easy to set as a digital clock.
- Medication is under lock and key, and is tamper resistant.
- Compartment divisions keep medications in place even if the unit is turned upside-down.
- MedReady uses both a battery and an A/C adapter during normal operation. The battery is used to dispense medication and the A/C power keeps the battery fully charged. When A/C power is removed, it will continue to dispense for 48 hours to allow for for travel or power failures and will charge back up when reconnected.
- "OFF" feature for long term disabling of the dispenser.
- "Early dosage" feature allows removal of medication prior to scheduled times.

MONITORED DEVICES


When the user does not respond to the medication alert in 30 minutes (adjustable), caregivers are notified immediately via email, text and phone calls. Full compliance reporting feature shows 2 years of medication history.

Cellular MR-357 Model

Compliance reporting in real time, change settings remotely, trigger early doses, change alarm times and more...

Landline 1750 Model

Compliance history sent once per day. Alarm times and other functions are set manually on the device.

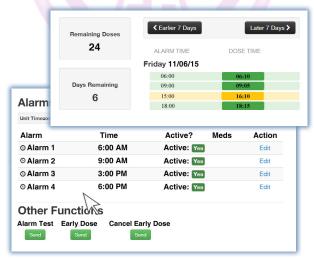

OPTIONS

Compliance Monitoring

- **1. Cellular Model MR-357.** Monitoring Fee. Includes all equipment and cellular plan. User does not need anything else.
- **2. Landline Model 1750.** Monitoring Fee. User needs non dedicated landline.
- 3. FL Alarm With
 Flashing Light
 for hearing impaired.
 Audible alarm PLUS
 Flashing Light.
- 4. 5310 Extra Pill Tray with clear plastic lid.
 Preload an additional tray for convenience.